

MOBIDERM®

A unique active substance in the treatment of œdema

- Foam blocks, separated by at least 1 cm, inserted between two dressings
- Difference in pressure exerted between the supporting zone and its surrounding area



Mobilisation of œdematous or indurated subcutaneous tissue following variations in muscle volume during activity

- Very rapid reduction of œdema



- Efficacy

- Speed of reduction in œdema volume
- Improved suppleness of fibrotic, sclerotic skin...



- Practicality

- Mobiderm, with its lightness, comfort and flexibility:
 - improves movements and hence drainage capacity
 - creates an efficient complementary treatment to drainage of the œdematous limb (does not dispense with the need for overall drainage, which stimulates the lymphatic system)
- The Mobiderm prefabricated system is :
 - easy to cut, making it easier for therapists to apply the bandage
 - washable and reusable (for up to 6 months by the same patient)



- Comfort

- extra-soft skin contact



Thuasne's Mobiderm range is a set of unique systems for the mobilization of oedematous or indurated subcutaneous tissues.



**Mobilizing
bandage**



**Mobilizing
pad**



Mobilizing garments

Sleeves for arm or leg
made-to-measure for
the patient



15 x 15 mm




5 x 5 mm

Dimensions:

	Mobilizing pad	Mobilizing bandage
Size	33 cm x 1 m	10 cm x 3m
blocks' size	5 x 5 mm 15 x 15 mm	15 x 15 mm

Principe of action :

Mobiderm creates a pressure differential between the supporting zone and the surrounding area, which, combined with muscle volume variations during exercise triggers a shear effect in the subcutaneous tissues concerned.

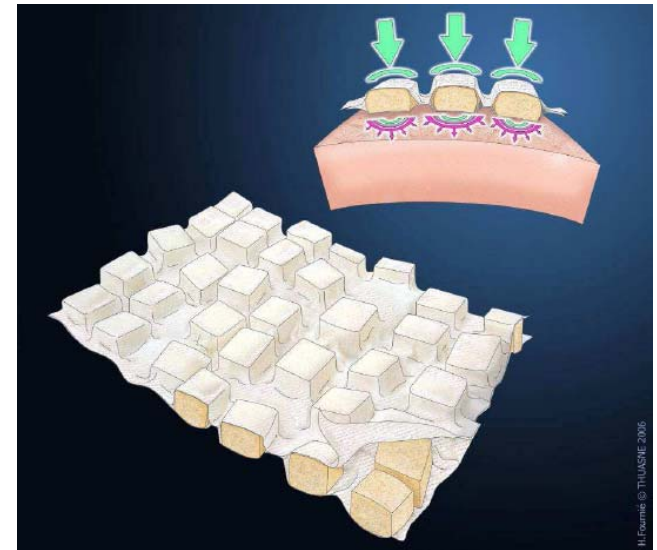

Composition:

adhesive non-woven fabric in polyester
Foam blocks in polyurethane

+

Polyamide fabric with 20% elasthanne



Bandage and pad



Made to measure
garment

Due to its action on the subcutaneous tissues and its shear effect Mobiderm® demonstrates encouraging results from the 5th day of treatment.



Mobiderm® treatment during 5 days

PRINCIPE ACTIF BREVETE POUR LE TRAITEMENT DES OEDEMES
(Lymphatiques / Veineux / Post-traumatique)



Before



3 months



After

ALFAcare



ALFAcare



Before

5 days

After

ALFAcare



Before

8 days

After

ALFAcare



Before



3 weeks

After

ALFAcare

MOBILIZING MADE TO MEASURE GARMENTS (Lymphoedema stabilisation phases)

Outside :

80% Polyamide
20% élasthanne

Inside:

56% Polyester (non tissé)
44% polyuréthane (cubes)



ALFAcare